

# SLEEPONE AUTO



SleepOne AutoCPAP device provides a comfortable therapy with the lowest possible therapy pressure. This is done by the internal algorithm. Reaction to Snoring is adjustable by the "Snore Sensitivity" Parameter. Small, lightweight and robust design makes it very portable together with universal power supply. Easy to adjust integrated heated humidifier, wide pressure range operation, auto on and off functions fulfill all the needs for a successful therapy.


Forget about costly accessories, you do not need to pay much for a decent therapy. No need to pay for getting your sleep data. Everything is included in this small, durable machine for a comfortable night sleep. A data logging function ensures you get the right therapy all the time. Bundled with the SleepOne Software package, you can get your reports with just connecting the device with a USB cable. Whether you use an AutoCPAP or a standard one; nightly reports could be fetched from the device with the aid of the SleepOne Software.

Internal power supply, ensures easy handling during travel and prevents the bulky mess in your bedroom. Only a thin power cable is enough to run the device. Optional accessories are available on demand. An oxygen enrichment adapter; a converter for use with 12V jack in the cars, trucks and boats;

SleepOne Platform gives the patient comfort, reliability and cost effectiveness.

## Technical Specifications

Pressure Range (cmH <sub>2</sub> O)	3 - 20 cm H <sub>2</sub> O (0.5 increments)
Flow	0 cm H <sub>2</sub> O 180 L/Min 20 cm H <sub>2</sub> O >50 L/Min
Altitude Adjustment	Automatic
Dimensions (WxDxH)	18.1 x 15.5 x 11.7 cm
Weight	1.28 Kg.
Power Consumption	65 Watt
Voltage	100-230 V AC
Warranty	2 Years
Noise Level	<29 dB
Report Software	SleepOne Software
Auto-On/Auto Off	Yes
Mask Off	Yes
Data management	SD Card stores up to 365 nights of detailed therapy data.
Adjustable Parameters	Ramp time, Ramp Starting Pressure, Auto-on, Auto-Off, CPAP Pressure low limit, CPAP Pressure high limit, Apnea parameter, Mask alert

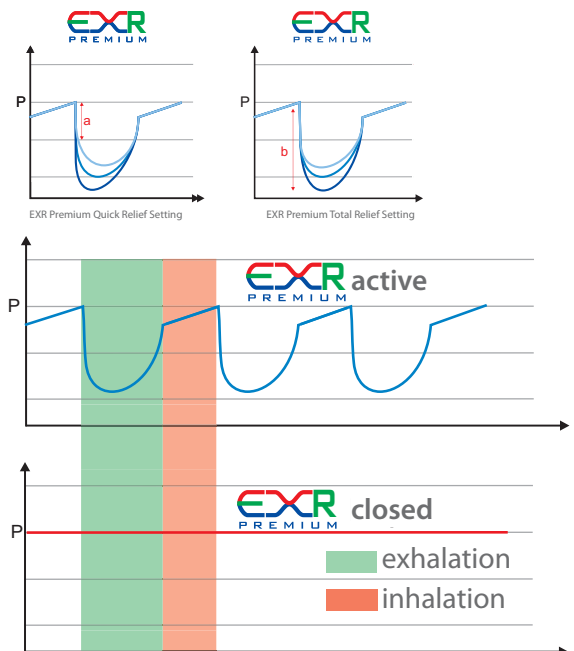

## Key Features

- Optimized algorithm enables high breathing comfort, allows ultimate patient compliance and increases the chance of therapy success
- Auto Leak Compensation
- Auto Altitude Compensation
- Low Noise Operation

